

STEPLITE® CL SHORT LEG WALKER

Model No: 43-316

Sizes: Small - Large

Color: Black



INDICATIONS:

- Post Cast Removal
- Severe Sprains
- Post Surgery
- Stable Fractures

BENEFITS:

Patients who need stability and protection during ambulation will benefit from the use of the StepLite®. This lightweight short leg walker contains a molded sole, side struts, a tri-laminate boot wrap, and five hook and loop straps providing a safe and secure fit. This device limits ankle motion and gives patients the ability to walk without the risk of further injury.

PRODUCT FEATURES:

- Made of a glass filled nylon material that is extremely lightweight while strong and durable
- Low heel height provides a smooth, energy conserving gait
- Boot wrap is made of C3™ tri-laminate that wicks moisture away from skin into Microban® treated foam for cool, clean and comfortable wear
- Easy and quick assembly

HOW TO APPLY:

Complete wearing instructions are enclosed in packaging.


SIZING:

Size according to shoe size.

Size	Men's Shoe	Women's Shoe
Small	6 - 7 1/2	5 - 8 1/2
Medium	8 - 9 1/2	9 - 10 1/2
Large	10 - 12 1/2	11 - 12 1/2

IMPORTANT:

When assembling the struts, they **must be snapped in far enough so that they are flush with the boot**. It might take a little extra strength to push the struts all the way in, but it is necessary to avoid any weak points.

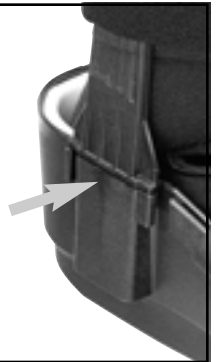

800-327-4110

FLA ORTHOPEDICS®

Website: www.flaorthopedics.com

©2002 FLA ORTHOPEDICS, INC.